

SINESON II

4 Macam Terapi dalam 1 Alat



SINESON II


4 macam terapi dalam **1** alat:

1. Extracorporeal Shockwave Therapy, jenis Pneumatik (ESWT)
2. Extracorporeal Shockwave Therapy, jenis Gelombang Piezo (ESWP)
3. Diatermi Frekuensi Radio (RF Diathermy): Resistif and Kapasitif
4. Diatermi Frekuensi Radio (RF Diathermy) Bipolar

PAIN RELIEF - PHYSIOTHERAPY - BODY CONTOURING

Extracorporeal Shock Wave Therapy, **Pneumatic Type (ESWT)**

Radial Type

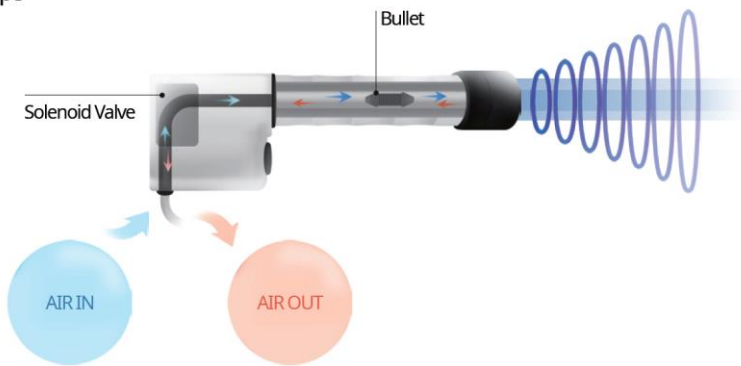


- Non-invasif
- Dalam tempo singkat mengurangi nyeri kronis, terutama yang disebabkan oleh deposit pengapuran (spur)
- Solusi untuk nyeri yang tidak respons terhadap terapi lain seperti anti-inflammatory atau immobilization devices



Extracorporeal Shock Wave Therapy, **Piezo Wave Type (ESWP)**

SINESON II

■ Applications | Urology and orthopedics


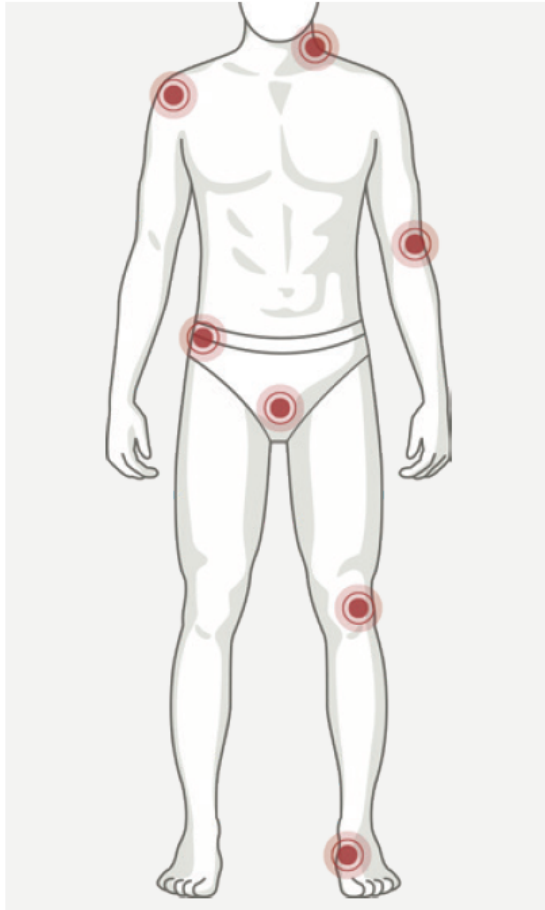


- Chronic neck, shoulder and back pain
- Shoulder calcifications, Rotator cuff tendonitis
- Chronic tenonitis, Chronic muscular pain
- Tennis or Golf elbow
- Myofascial pain syndrome, Trigger point syndrome

- Patellar tendonitis, Knee osteoarthritis
- Achilles tendon pain
- Heel pain, Plantar fasciitis
- Degenerative arthritis, Rheumarthritits

- Prostatodynia
- Prostatitis (All kinds of infection related to prostate gland)
- Prostate stone
- Prostate enlargement
- Inflammation disease of Anal fistula or Perianal abscess



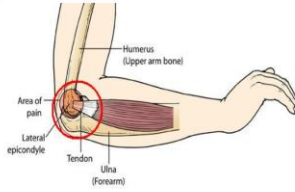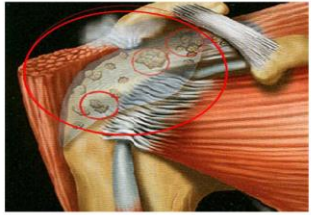

Treat 3~5 times on areas marked with a red circle in above figures for 3~5 Min.

●●● Inflammation Pain Care with ESWP

⌚ 2 sessions per week Total 6 sessions / 1 session – 5 min.





●●● Prostatitis, Prostate stone Care with ESWP

⌚ 2 ~ 3 sessions per week total 20 sessions / 1 session – 5 min.
Like lipocyte, prostate be reduced by shock.

[Prostatitis]

Care with ESWP

[Normal prostate]

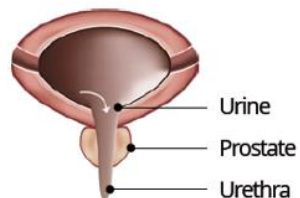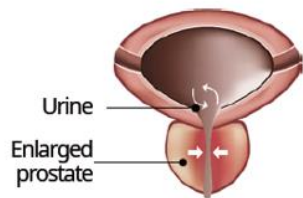

■ Postures for BPH treatments

· 01





· 02




· 03



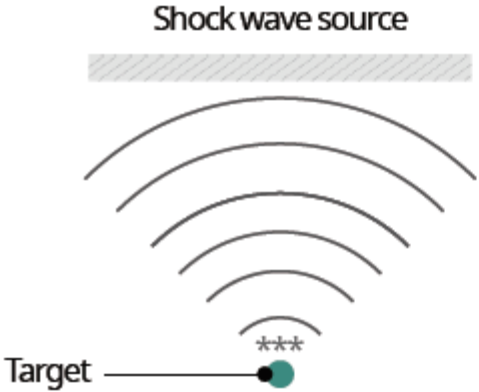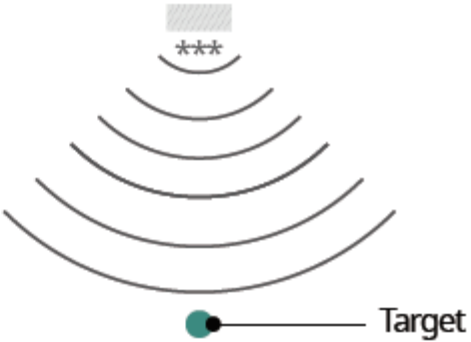

ESWP	ESWT
Focused extracorporeal shock waves	Radial extracorporeal shock waves
Shock wave source  Target	Shock wave source  Target
*** Point of maximum energy flux dencity	

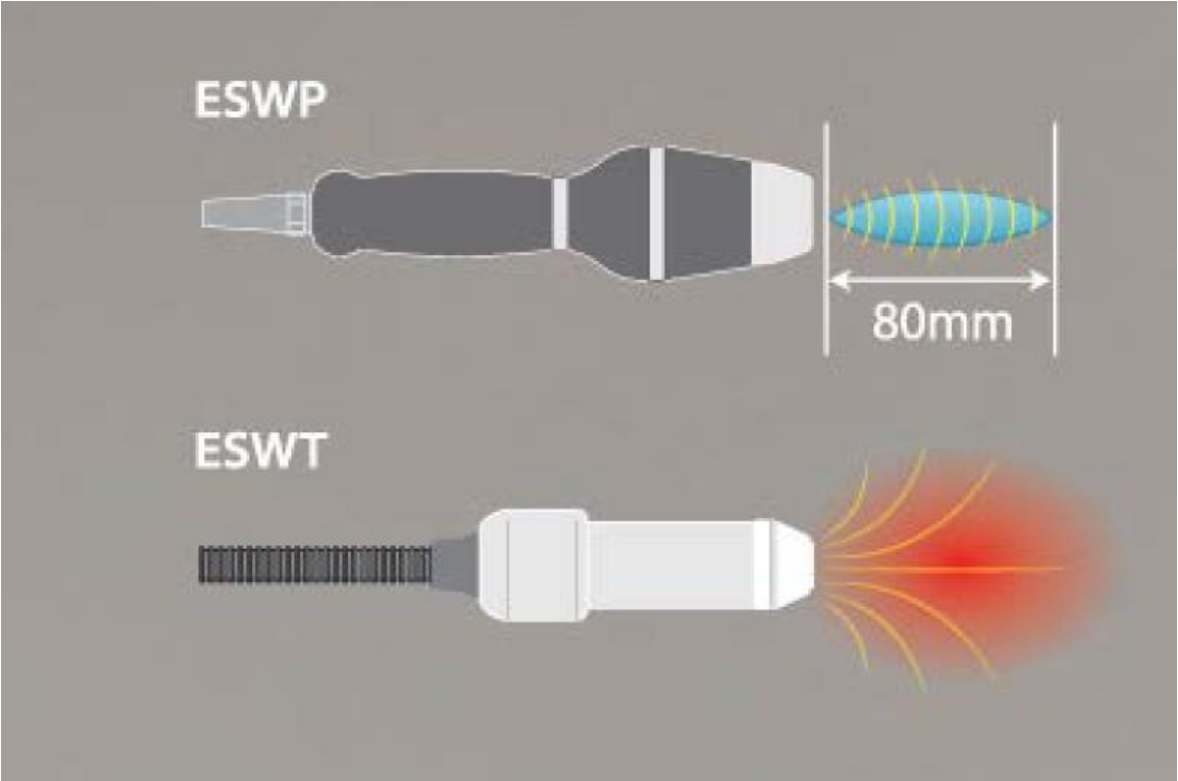




RF (Radio Frequency) Diathermy

Accessories:

- RET electrodes (Resistive Electrode Transfer)
- CET electrodes (Capasitive Electrode Transfer)
- Bipolar electrodes
- Sineston set


Accessories:

- RET electrodes 
- CET electrodes 
- Bipolar electrodes
- Sineston set

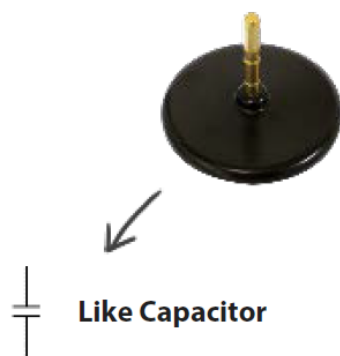

For Soft tissue care

CET Capacitive Electrode Transfer

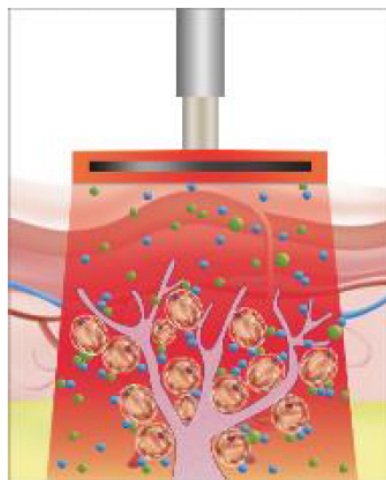


Face care, Hair loss care,
Subcutaneous fat care, mussels etc.

Non conductive electrode – Insulated electrode
INDUCED CURRENT

Hard tissue care

RET Resistive Electrode Transfer

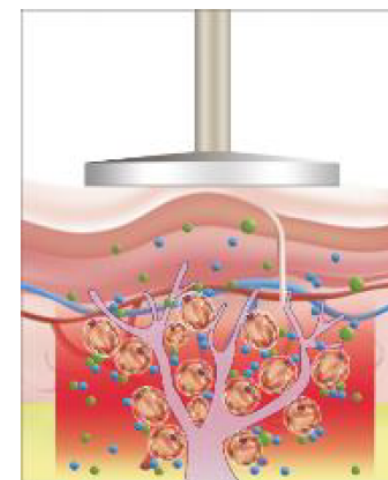


Joint, tendons, bones, Back,
Visceral fat,
Cellulite care etc.

Conductive electrode – Non insulated electrode
Direct CURRENT with less dispersion of energy
Go into deeper than CET but narrow part


Accessories:

- RET electrodes 
- CET electrodes 
- Bipolar electrodes
- Sineston set

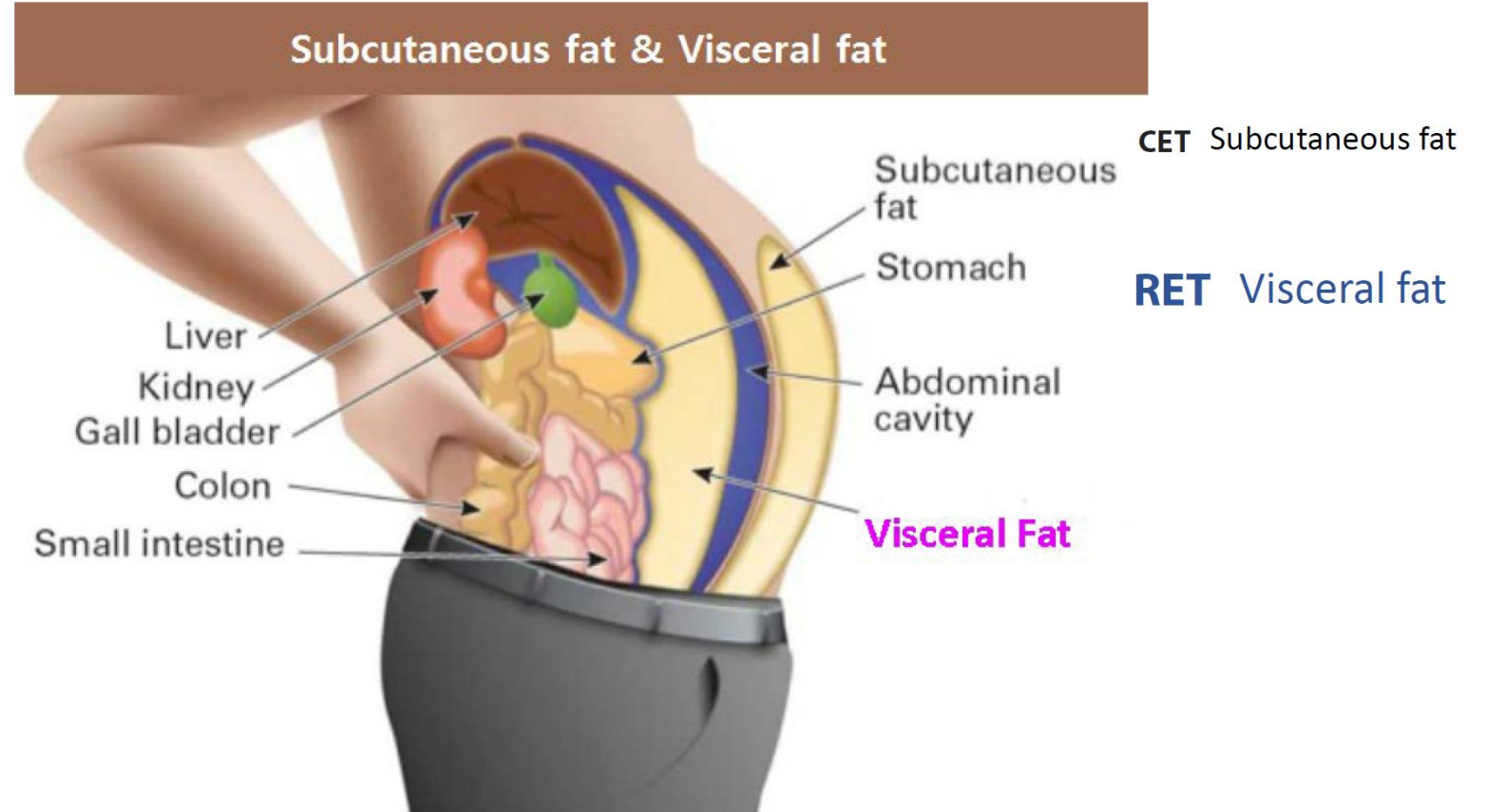

Accessories:


- RET electrodes 
- CET electrodes 
- Bipolar electrodes
- Sineston set

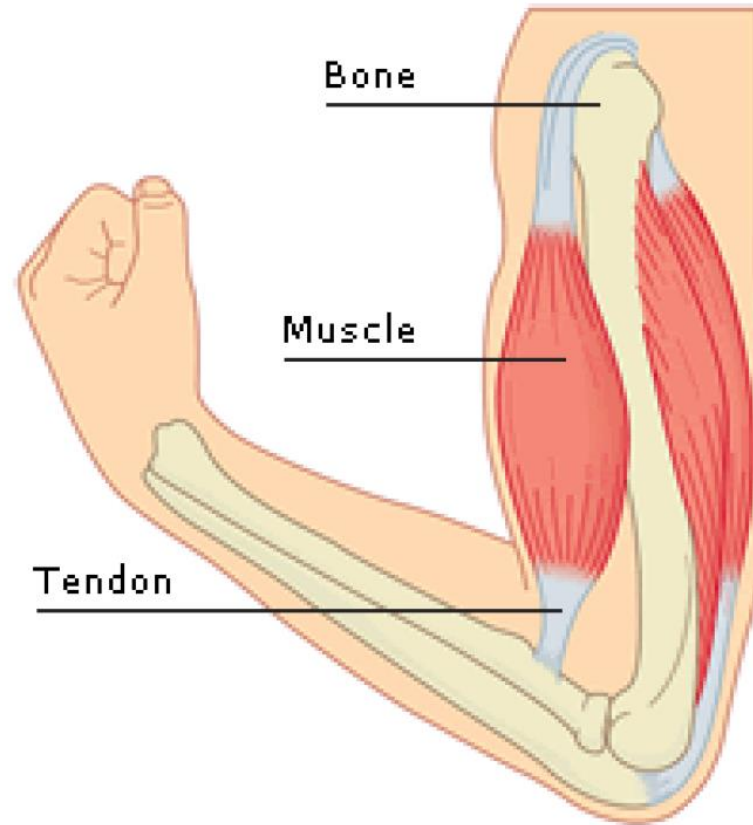

- CET** For Soft tissue care
-Subcutaneous fat, Muscles
- RET** Hard tissue care
- Joint, tendons, bones

Accessories:


- RET electrodes
- CET electrodes
- Bipolar electrodes
- Sineston set

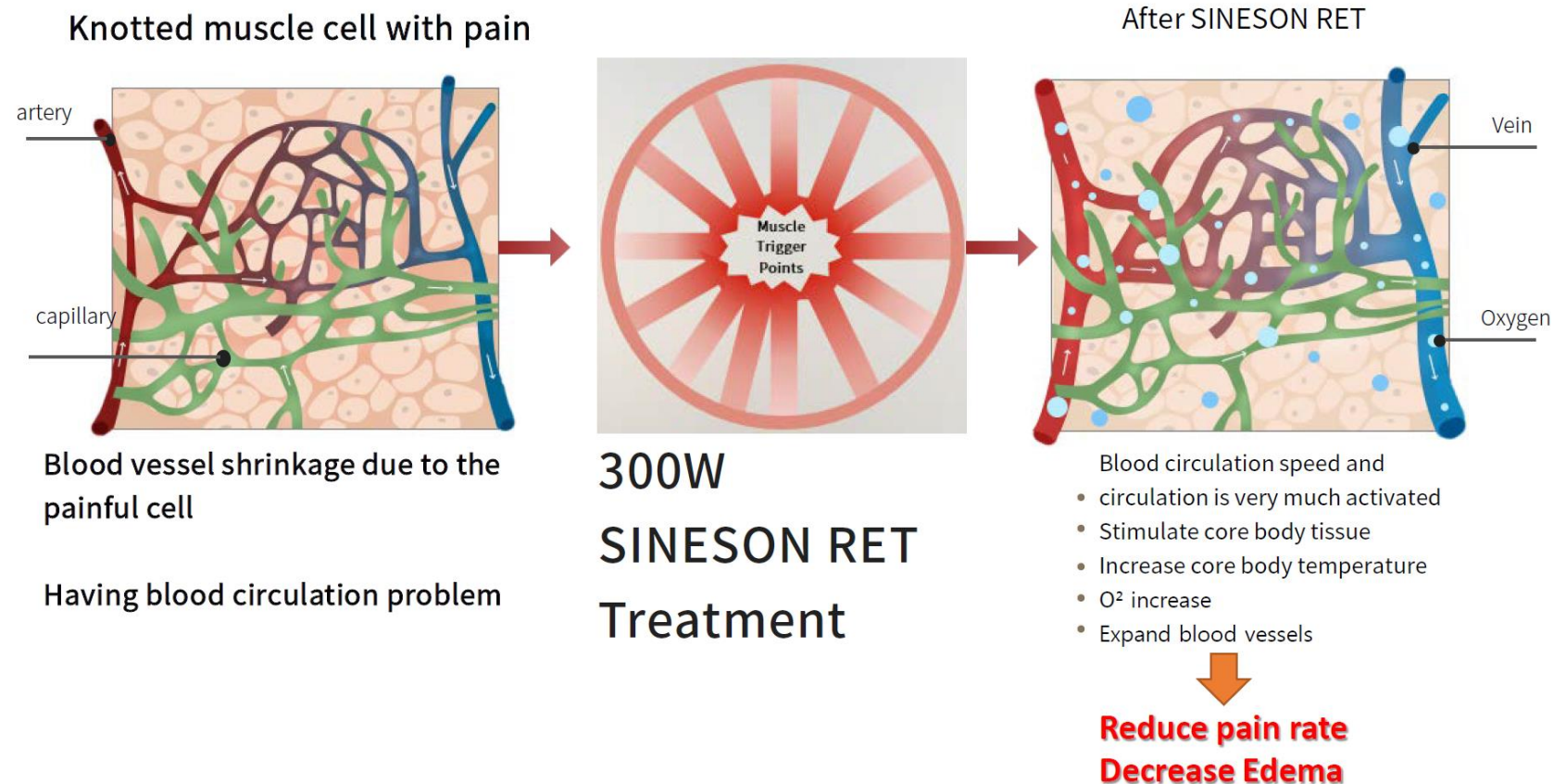

Accessories:

- RET electrodes →
- CET electrodes
- Bipolar electrodes
- Sineston set

EFFECT ON EYE

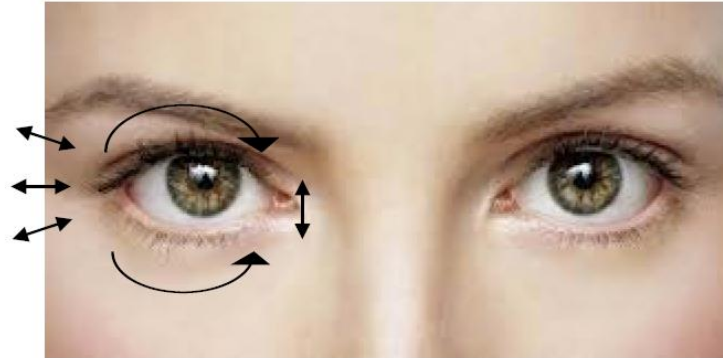

- Dark circle
Expediting generation of cells by enough blood inflow to eyelids and lymph circulation, eye-care discharges toxic matters and removes the dark circle.
- Eye edema
Eye edema by accumulation of fatigue and tired eyes by eye muscle tension and hypofunction of lacrimal glands can be cured by the eye-care.

RF (Radio Frequency) Diathermy


SINESON II

Accessories:




- RET electrodes
- CET electrodes
- **Bipolar electrodes**
- Sineston set



Mono Polar RF

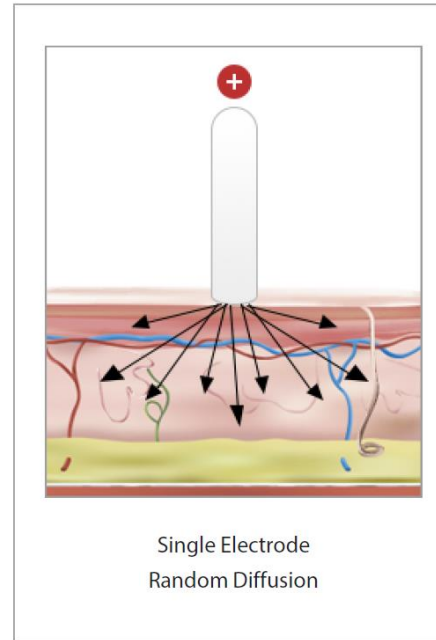


Bi Polar RF

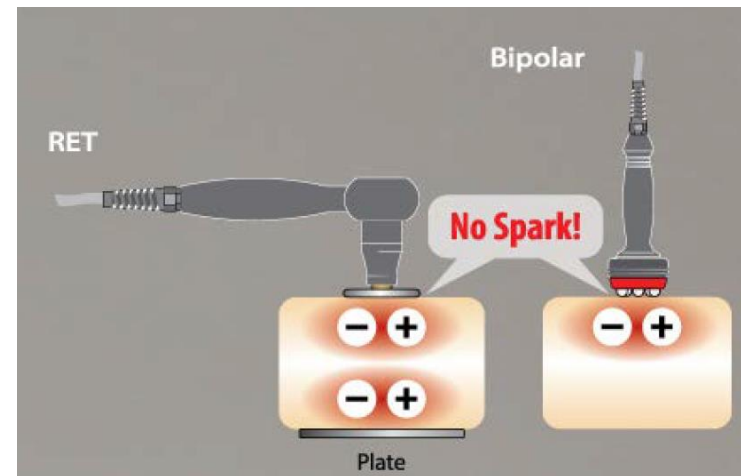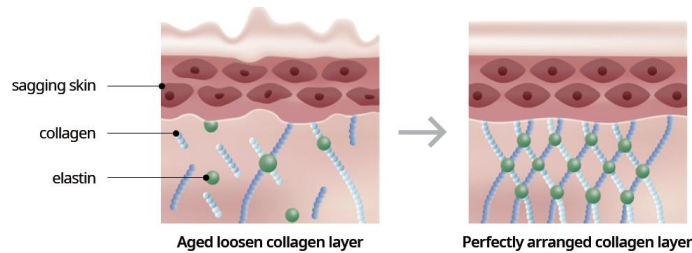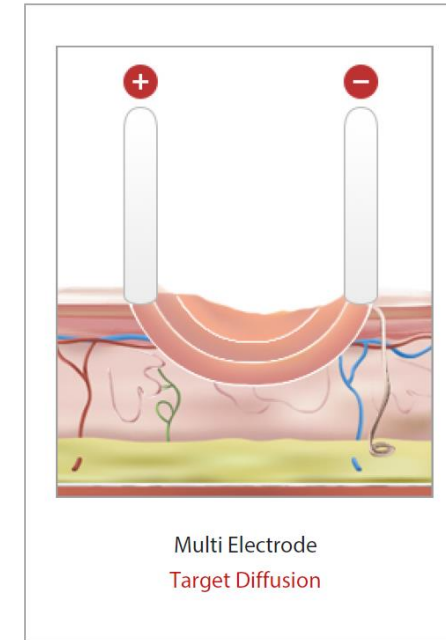

Accessories:

- RET electrodes
- CET electrodes
- Bipolar electrodes
- **Sineston set** 

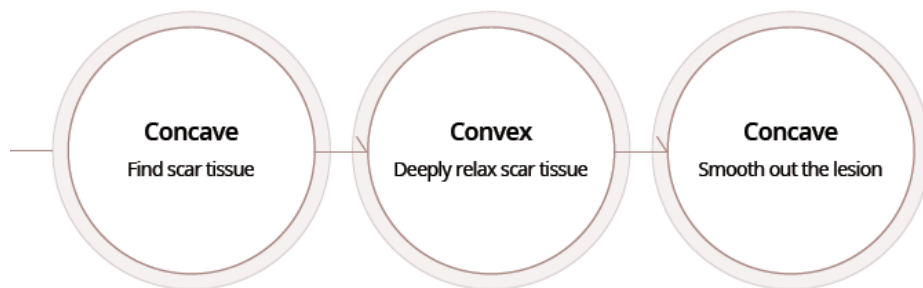

- | | |
|--|---|
|  Concave | - Equalizes or disperses pressure over a larger surface area
- Useful for increasing comfort, especially with a pain-sensitive patient |
|  Convex | - Maximize pressure in a small surface area
- Allows clinician to pinpoint pressure accurately |



Accessories:

- RET electrodes
- CET electrodes
- Bipolar electrodes
- **Sineston set** 

MFR(Myofascial Release with Diathermy Transfer)



RELAX FASCIA, REDUCE PAIN WITH SINESTON - OPTIONAL ITEM




The benefits of SINESTON

- 01 Increase core body temperature with RF
- 02 Comfortable feeling with diathermy without unstable time at the moment of first contact with tool on the body.
- 03 Faster effect thanks to synergy effect with RF than tool use only
- 04 Shorten treatment time
- 05 Less fatigue than meridian scrapping massage




RF (Radio Frequency) Diathermy

SINESON II

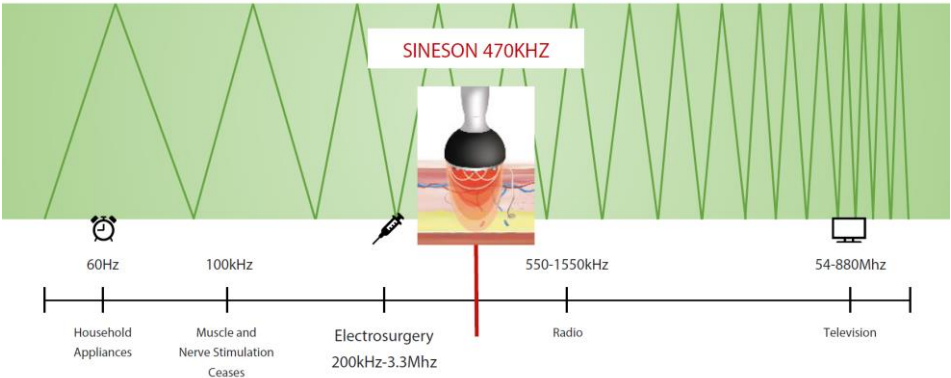


EFFECT ON SKIN

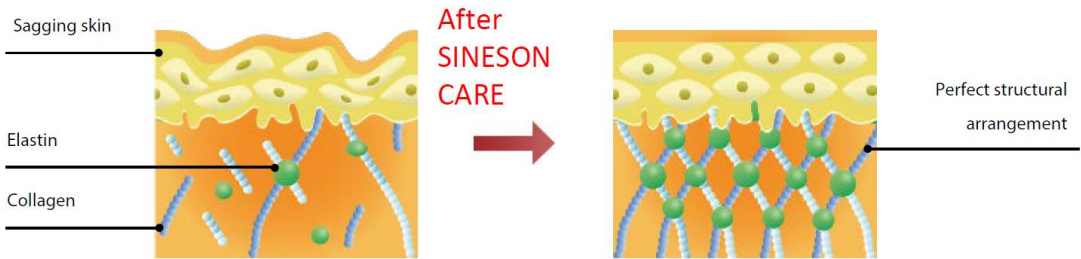

Damage old weak collagen
↓
Activate a generating neo collagen and elastin
↓
Skin tightening, lifting, thicker dermis layer effect

EFFECT ON CELLULITE

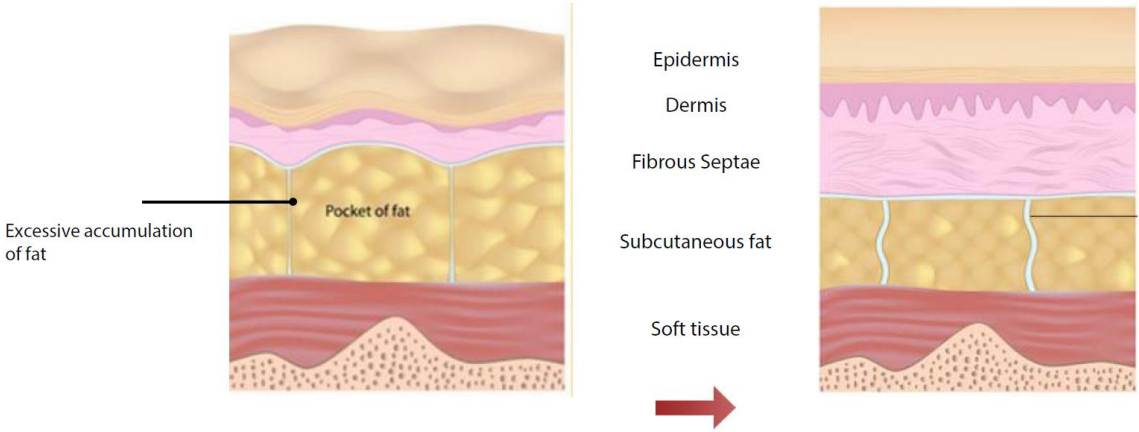

SINESON II



- 1. 300W High Power RF** enables all types of clients and all kinds of face, body and physiotherapy treatment and guaranty the maximize RF result
- 2. Technology intensive Power System to minimize power loss**
 - Installed P.F.C
 - LLC Current Resonant Converter
- 3. High durability installing protection circuit**
 - O.C.P._Over Current Protection
 - O.V.P. Over Voltage Protection
 - O.L.P._Over Load Protection
 - A.O.C.P._Abnormal action Over Current Protection
 - T.S.D _Internal Thermal Shutdown
 - O.S.S.P._Output Short Shutdown Protection
- 4. Comfortable SPARK FREE treatment installing spark free detector**
- 5. User friendly user interface**